

prosthetic.line

## CM LOC<sup>®</sup> Post coping anchor

|    |                      |          |    |
|----|----------------------|----------|----|
| DE | Gebrauchsanweisung   | Deutsch  | 1  |
| FR | Mode d'emploi        | Français | 12 |
| EN | Instructions for Use | English  | 22 |
| IT | Modo d'uso           | Italiano | 32 |
| ES | Instrucciones de uso | Español  | 42 |
| DA | Brugsanvisning       | Dansk    | 52 |
| ZH | 使用说明书                | 简体中文     | 62 |
| ZH | 使用說明書                | 繁體中文     | 72 |
| JA | 取扱説明書                | 日本語      | 82 |
| KO | 사용 설명서               | 한국어      | 92 |

# Instructions for Use CM LOC® Post coping anchor

## 1 Scope of application of Instructions for Use

These Instructions for Use apply to the products listed under Section 29. The issuing of these Instructions for Use renders all previous versions invalid. The manufacturer rejects any liability for damages resulting from non-compliance with these Instructions for Use.

## 2 Trade name

See Section 29.

## 3 Intended use




The products are intended for prosthetic restorations and to support procedures in the dental clinic or laboratory.

## 4 Expected clinical benefit

Restoration of chewing function and improved aesthetics.

The Summary of Safety and Clinical Performance, SSCP for the implantable devices covered by these Instructions for Use, is available on our website and accessible at this address: [www.cmsa.ch/docs](http://www.cmsa.ch/docs).

## 5 Product description

| Product  | Description  |
|--|--|
|   | <b>CM LOC® Male parts</b><br>Castable, solderable or laserable male parts for hybrid-supported, removable dental prostheses on root canal caps.                    |
|   | <b>CM LOC® Female part</b><br>Retaining element as connecting part between denture and abutment.   |
|  <p>                         extra-low<br/>ca. 400g                              low<br/>ca. 1200g                              medium<br/>ca. 1800g                              strong<br/>ca. 2400g                     </p> | <b>CM LOC® Retention inserts</b><br>Exchangeable retention inserts in four defined force levels.<br>yellow: extra-low<br>red: low<br>green: medium<br>blue: strong |

## 6 Indications

### CM LOC® male part C and CM LOC® male part E

Implant anchorage for hybrid-supported removable dental prostheses on root canal caps, in combination with the specific CM LOC® system for female parts.

#### Mandible

Anchorage of mandible (MD) prosthesis on two or more root canal caps.

#### Maxilla

Anchorage of maxillary (MX) prosthesis on four or more root canal caps.

## 7 Contraindications

- Divergences > 20° (per anchoring element).
- Use on a single root canal cap.
- Restoration of severely periodontally damaged abutment teeth.
- Not suitable if a fixed connection is required.
- Unilateral free-end prosthesis without transversal support.
- Patients who are unable to keep the regularly required check-up appointments for health reasons.
- Patients with bruxism or other para-functional habits.
- Patients with allergies to materials used in the product, see Section 19.
- Existing clinical picture in the patient's mouth does not permit the correct application of the products.

## 8 Compatible products

The CM LOC® system for female parts is also compatible with the following Locator®-like male parts:

- CM LOC®
- CM LOC® FLEX
- MedentiLOC®

Please contact us for further information regarding compatibilities.




The retention force on these abutments may vary due to the different manufacturing tolerances and surfaces of the various male parts.

To fabricate the finished denture, a number of general laboratory supplies are required in addition to the products listed under Section 29. The following gives a selection of materials that Cendres+Métaux SA offers in its portfolio.

|          |                         |          |                                     |
|----------|-------------------------|----------|-------------------------------------|
| 08052138 | Polyurock Kit           | 08055014 | Livento® invest Powder (50 x 100 g) |
| 08052135 | Polyurock Catalyst      | 083739   | Livento® invest Liquid (1000 ml)    |
| 08052136 | Polyurock Release Spray | 08052160 | uniVest® Plus Powder (30 x 150 g)   |
| 08052137 | Polyurock Mixer         | 08052161 | uniVest® Plus Liquid (1000 ml)      |
| 08052566 | Polyurock Colour yellow | 08052162 | uniVest® Rapid Powder (30 x 150 g)  |
| 08052149 | ABF Wax Universal       | 08052163 | uniVest® Rapid Liquid (1000 ml)     |
| 08052150 | ABF Wax Creativ light   | 080181   | CM soldering investment (4 kg)      |
| 08052151 | ABF Wax Creativ dark    | 080229   | CM soldering paste                  |
| 08052154 | ABF Wax Special         | 08052307 | Legabril Diamond (50 g)             |
| 08052148 | ABF Wax Margin          |          |                                     |
| 08052153 | ABF Wax Position        |          |                                     |
| 08052152 | ABF Wax Tecno           |          |                                     |

## 9 Qualification of the specialist

Expertise in professional dentistry and dental technology is assumed. The current Instructions for Use must be available at all times and be completely read and understood before the first application. The fabrication of dentures and their maintenance may only be performed by qualified specialists.


 Important information for the specialist

 Warning symbol for increased caution

## 10 Prescription

Federal laws in the USA prohibit the use by or sale to unlicensed dentists.

## 11 Side effects

 This product must not be used in patients with allergies or suspected allergies to materials used in the product (see Section 19), or only after prior allergological clarification.

Auxiliary instruments may contain nickel.

If applied as intended, side effects can be excluded.

## 12 Warnings

### Magnetic resonance (MR) environment

The device has not been evaluated for safety and compatibility in the MR environment.

The product has not been tested for heating or migration in the MR environment.

### CM LOC® Spacer

The CM LOC® Spacer is slightly oversized with regard to the original components.

This ensures optimal space conditions for later polymerisation in the mouth.

The spacer must not be used in place of the female part or as a temporary replacement.

## 13 General information



– Wax-up of the root canal cap with root canal post: If there are multiple root canal caps, prepare the solder/laser surface at right angles to the direction of insertion. Use prefabricated, cast-on precious metal pins.

– After soldering/casting, slowly cool to room temperature. The optimal mechanical properties are achieved allowing it to bench-cool to room temperature. Fit the CM LOC® spacer to protect the male part while sandblasting and processing.

– The height of the CM LOC® Block-out spacer can be shortened to achieve better mounting of the female part.

– When working with the CM LOC® male part for retention of overdentures, the direct or indirect approach may be used.

– We recommend that the clinical case be designed such that the largest possible support polygon is achieved. Small distances between consecutive implants and long free-end saddles can cause undesirable effects such as increased wear of the system components.

– Proper seating of the dentures on the mucosa must be checked at least once each year, and relining carried out if required to prevent rocking movement (overload). We recommend checking the denture at regular intervals of approx. 3 months and to replace the retention inserts if necessary.

– When fabricating new dentures and when using a palate-free design, we recommend fabricating an individual reinforcement framework.



– In patients with suspected titanium allergy or hypersensitivity, we alternatively recommend the use of the Pekkton® female part. One must allow for an increased aftercare effort and, if necessary, changing/replacing the system for female parts, as Pekkton® is somewhat softer than titanium as the material for the female part.

### Integration of the housing of the female part


#### Direct method

The dentist providing treatment may integrate the CM LOC® Housing of the female part and the retention inserts in an existing or new prosthesis directly during the treatment session.

#### Indirect method


The dentist must take an impression of the CM LOC® male part with the CM LOC® Impression Part and send the impression to the laboratory for subsequent fabrication of the model. The laboratory then inserts the CM LOC® Analog in the CM LOC® Impression part to facilitate reliable transfer of the position of the CM LOC® Male part in the mouth, and fabricates the master model.

**14 Preventive measures**


-  – The product components are supplied non-sterile. For more information see Section 16 "Reprocessing".
- Only original tools and parts may be used for this work. For information and additional details, please contact your Cendres+Métaux SA representative.
- Before any procedure, ensure that all required product components are available in sufficient quantity.
- For your own safety, always wear suitable protective clothing. In particular when grinding, we recommend wearing protective goggles and a dust mask as well as the use of a suction unit.
- Secure parts against aspiration.
- The mechanical cleaning by patients with a toothbrush and toothpaste may lead to premature wear.

**15 Single use**

Products that are intended for single use and are labelled "single-use" accordingly are subject to a certain amount of stress, increased wear, and even loss of functionality during their use.

-  Multiple application of products labelled «single use» was not tested. This can impair the safety, function and performance of the products as well as increase the risk of transmitting infections.

**16 Reprocessing**

-  The prosthetic work, including all system components, must be cleaned, disinfected and, if appropriate, sterilised prior to each work step. Materials made of metal alloys, high-performance polymers (Pekkton®) and ceramics are suitable for steam sterilisation. With the exception of Pekkton®, components made of plastics are not suitable for steam sterilisation. Consider published national guidelines when selecting a disinfection and sterilisation process and the Instructions for Use "Reprocessing of surgical and prosthetic products" ([www.cmsa.ch/docs](http://www.cmsa.ch/docs)).

**17 Scope of application**

After fabricating the root post cap, the CM LOC® male part can be inserted or processed according to the indication. The CM LOC® male parts are designed to fixate overdentures (total dentures) or partial dentures completely or partially through root canal caps in the maxilla and/or mandible using the specific CM LOC® System for female parts. The CM LOC® Anchoring system consists of a root canal cap consisting of two standardised male parts for casting, soldering or lasering with which implant divergences of up to 20° are compensated and a system for female parts with four exchangeable retention inserts in four defined force levels.

**18 Procedure**


**Inserting the CM LOC® male part E in Elitor® by laser welding.**

Initial situation.




**Preparation for laser welding**

In a first step, face mill the root canal cap already fabricated using a milling machine at right angles parallel to the direction of insertion.

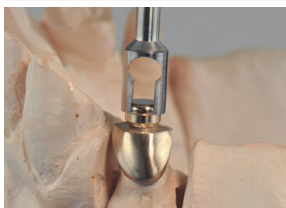


With the CM LOC® Parallelometer insert, set the CM LOC® male part as centrally as possible on the already cast, face-milled root canal cap and securely fix it with wax using the root canal cap.

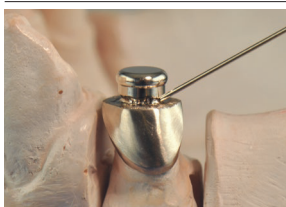

**Laser welding**



Then fill all undercuts circularly around the entire CM LOC® male part in the laser unit with laser wire.

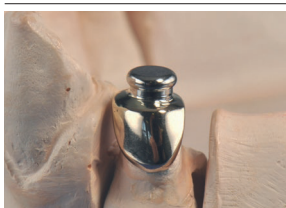-  Please observe the manufacturer's instructions on the laser procedure.



**Finishing**

Then, using a standard rubber wheel, smooth down the laser welds and polish using a polishing brush.

-  To simplify working with and protecting the CM LOC® male part, put the CM LOC® Spacer on the CM LOC® male part.
-  Make sure that no more material is removed than to the outer bottom edge of the CM LOC® male part.

**Inserting the CM LOC® male part in Ceramicor® by casting-on.**

Initial situation.



**Preparation for casting-on**


With the CM LOC® Parallelometer insert, set the CM LOC® male part as centrally as possible on the already modelled root canal cap and securely fix it with wax using the root canal cap.




**Casting / Finishing**


Then embed and cast.

Please observe the instructions for use for the dental casting alloys [www.cmsa.ch/docs](http://www.cmsa.ch/docs)

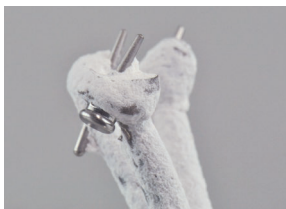


Fit the CM LOC® Spacer to protect the CM LOC® male part while sandblasting and processing. Then, using a standard rubber wheel, smooth down the cast and polish using a polishing brush.

 To simplify working with and protecting the CM LOC® male part, put the CM LOC® Spacer on the CM LOC® male part.

 Make sure that no more material is removed than to the outer bottom edge of the CM LOC® male part.

**Inserting the CM LOC® male part in Ceramicor® by soldering-on.**

Initial situation.

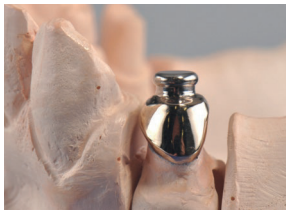

**Preparation for soldering-on**

With the CM LOC® Parallelometer insert, set the CM LOC® male part as centrally as possible on the already cast, face-milled root canal cap and securely fix it circularly with wax using the root canal cap.




**Soldering**

The solder gap should be continuous and between 0.05 – 0.20 mm wide.


Then design the soldering block, so that the male part is securely held and good flame access is ensured (observe corresponding solidus).




 – After soldering, slowly cool to room temperature. The optimal mechanical properties are achieved allowing it to bench-cool to room temperature.

**Finishing**

Then, using a standard rubber wheel, smooth down the solder joint and polish using a polishing brush.

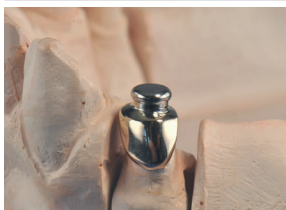

 To simplify working with and protecting the CM LOC® male part, put the CM LOC® Spacer on the CM LOC® male part.

 Make sure that no more material is removed than to the outer bottom edge of the CM LOC® male part.

**Impression taking of the oral situation for further processing using the indirect method**




Place the CM LOC® Impression part on the CM LOC® male part and create a functional impression. Ensure that the CM LOC® Impression part is correctly seated. Use a solid impression material (e.g. Impregum™).



**i** Check that the material is fully distributed around the CM LOC® Impression part and that no impression material has spilled into the CM LOC® Impression part.

Otherwise, clean the abutment and repeat the impression-taking process.

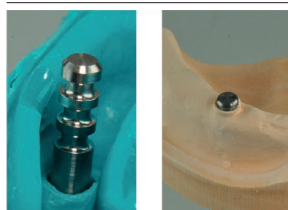

Then pass to dental laboratory for fabrication of the model. To fabricate the model in the laboratory, place the CM LOC® Analog in the CM LOC® Impression part and fabricate the master model.



Then place either the CM LOC® Housing of the female part with a mounted CM LOC® Processing insert or the CM LOC® Spacer onto the CM LOC® Analog. Use of the CM LOC® Spacer or the original CM LOC® Housing of the female part is at the discretion of the user.

**i** Spacer = place holder for polymerisation in mouth.  
Direct polymerisation with CM LOC® Housing of the female part in the laboratory.  
Block out all undercuts during further processing and fabrication of the denture.




**Finalisation**

The prosthesis can now be fabricated using conventional technology. After processing, remove excess resin around the CM LOC® Housing of the female part with a round bur. Then finish and polish the prosthesis. Afterwards replace the processing insert in the CM LOC® Housing of the female part with a Pekkton® Retention insert in the desired force level.

**i** See description in Selection of retention inserts.

**i** Make sure that no polymer has flowed into the housing of the female part. If necessary, remove the processing insert and carefully remove excess polymer from the inside of the housing of the female part with a probe.

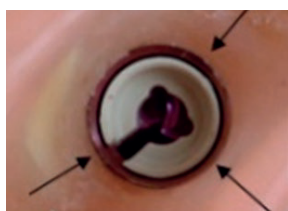

**Selection of retention inserts**


Four different CM LOC® Retention inserts made of Pekkton® are available for retention. The retention inserts are colour-coded and divided into four different levels of retention force.

yellow: extra-low    red: low  
green: medium    blue: strong

**i** Ensure that the selection of pull-off forces is adapted to the clinical situation. Only use the extra-low insert for immediate loading at the start.

**i** To provide the patient with comfortable and easy integration of the prosthesis as well as familiarisation with the retention in the mouth, it is recommended to fit the prosthesis with a CM LOC® Retention insert, extra-low, first. If the patient demands stronger retention, CM LOC® Retention inserts offering greater retention force can be used. For assembly and disassembly of the retention inserts, see description: assembly and disassembly of retention inserts.

**Direct method: processing of the CM LOC® Housing during the treatment session.**

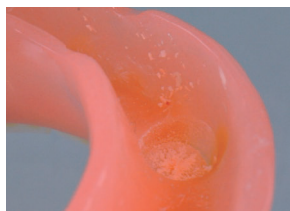

It is essential to create sufficient space in the prosthesis prior to inclusion in the prosthetic body. Use a standard round bur for this purpose. There must not be any contact between the prosthesis and the CM LOC® Housing of the female part.




Mount the CM LOC® Block-out spacer on the male part.

**i** Make sure that the CM LOC® Block-out spacer fits well. The height of the CM LOC® Block-out spacer can be shortened to achieve better mounting of the female part.




Then mount the CM LOC® Housing of the female part with mounted processing insert on the male part.

**i** Make sure that all undercuts are blocked out before polymerisation. Use a cold-curing polymer (e.g. GC Reline™, GC Advanced Technologies® Inc.) to anchor the CM LOC® Housing of the female part in the prosthesis. Apply the cold-curing polymer in the exposed area in the prosthesis and around the CM LOC® Housing of the female part.

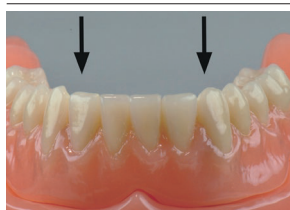

**Finalisation**

Place the prosthesis on the CM LOC® Male part in the oral cavity. Make sure that the prosthesis is entirely in occlusion with the opposing jaw. Ensure that the prosthesis is retained passively without compression on the soft tissue while the cold-curing polymer cures. Excessive occlusal pressure during curing can cause the soft tissue to be compressed and then decompressed again. This can cause the processing inserts to then click out of position.

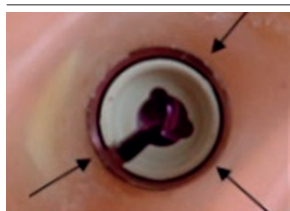

After processing, take the CM LOC® Block-out spacer out of the mouth. After processing, remove excess resin around the CM LOC® Housing of the female part with a round bur. Then finish and polish the prosthesis.

Afterwards replace the processing insert in the CM LOC® Housing of the female part with a Pekkton® Retention insert in the desired force level.

**i** See description in Selection of retention inserts.

**i** Make sure that no polymer has flowed into the housing of the female part. If necessary, remove the processing insert and carefully remove excess polymer from the inside of the housing of the female part with a probe.



**Selection of retention inserts**

Four different CM LOC® Retention inserts made of Pekkton® are available for retention. The retention inserts are colour-coded and divided into four different levels of retention force.

yellow: extra-low    red: low  
green: medium    blue: strong

**i** Ensure that the selection of pull-off forces is adapted to the clinical situation. Only use the extra-low insert for immediate loading at the start.


**i** To provide the patient with comfortable and easy integration of the prosthesis as well as familiarisation with the retention in the mouth, it is recommended to fit the prosthesis with a CM LOC® Retention insert, extra-low, first. If the patient demands stronger retention, CM LOC® Retention inserts offering greater retention force can be used. For assembly and disassembly of the retention inserts, see description: assembly and disassembly of retention inserts.

**Assembly and disassembly of the retention inserts.**




**Assembly**


The retention inserts are placed in the housing of the female part using the CM LOC® Multi-tool provided. Pick up the CM LOC® Retention insert with the IN side.

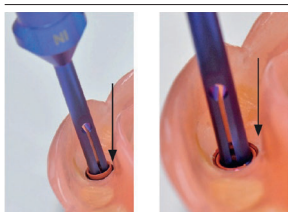

You can feel and hear the CM LOC® Retention insert lock into place. Press the CM LOC® Retention insert into the CM LOC® Housing of the female part in straight and parallel fashion until you can feel and hear it click into place.

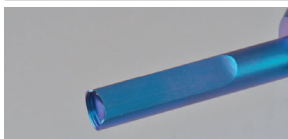


**Disassembly**  
Position the CM LOC® Multi-tool with the OUT side straight and parallel in the space between the CM LOC® Housing of the female part and the retention insert and press lightly into the CM LOC® Housing of the female part. The CM LOC® Retention insert thus unlocks and can be removed in a straight manner from the CM LOC® Housing of the female part.

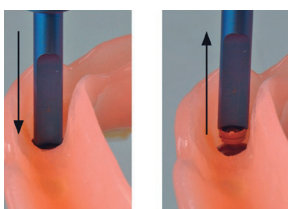


**Disassembly CM LOC® Housing of the female part.**



To disassemble the CM LOC® Housing of the female part, use the CM LOC® Extractor for the housing of the female part.

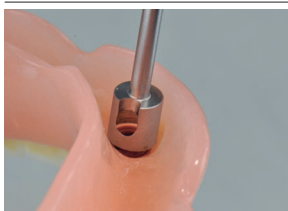

Mill the complete CM LOC® Housing of the female part with the CM LOC® Extractor for the housing of the female part.

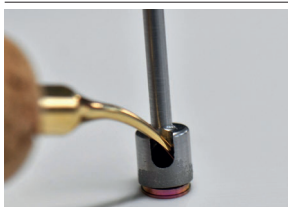

Then remove the CM LOC® housing of the female part from the CM LOC® Extractor for the housing of the female part through the side opening with an instrument. For better withdrawal, it is recommended to briefly warm the CM LOC® Extractor for the housing of the female part over a flame.

**19 Materials**

**C = Ceramicor®;** Au 60.0%, Pt 19.0%, Pd 20.0%, Ir 1.0%.

$T_s - T_L$  1400 – 1490 °C.

**E = Elitor®;** Au 68.6%, Pt 2.4%, Pd 3.9%, Ag 11.8%, Cu 10.6%, Zn 2.5%.

$T_s - T_L$  880 – 940 °C.

More detailed information on the materials as well as their compositions can be found in the product-specific material data sheets, the product information as well as the product list compiled in Section 29. All relevant documents can be found on the website [www.cmsa.ch/docs](http://www.cmsa.ch/docs) by entering the relevant product name.

**20 Notes on storage**



Insofar as no specific information on storage is given on the packaging of the product, we recommend storing the product in its original packaging, in a dry place, at room temperature and without direct sunlight. Improper storage can influence the product properties and lead to failure of the restoration.

**21 Patient information**

**21.1 Handling / follow-up**

On the day of insertion of the dentures at the latest, the patient must be informed that regular follow-up care is necessary to maintain the health of the entire masticatory system and the functionality of the denture. Ensure that the patients are motivated and instructed with regard to caring for their teeth as well as dentures.

Permanent and removable dentures are subject to considerable stress. Signs of wear are normal and cannot be avoided, only reduced. The amount of wear depends on the overall system.



Our endeavours are aimed at using materials that are as optimally matched as possible in order to reduce wear to an absolute minimum. Proper seating of the dentures on the mucosa must be checked at least once each year, and relining must be performed if required to prevent rocking movement (overload). We recommend checking the dentures at intervals of approx. 3 months initially and to replace the auxiliary parts such as retention inserts if necessary.

## 21.2 Insertion and removal of the dentures

It should be ensured that the dentures do not tilt, as any tilting can lead to damage. The denture should never be inserted by clenching the teeth, as this can damage or even break the connecting element.

### Insertion

The denture can be placed on the anchor elements in the mouth using the thumb and index finger. Then it is correctly positioned on the anchoring elements applying gentle, even pressure. By carefully closing the jaws, it is possible to check whether the denture is in its correct final position.

### Removal

For removal, the denture can be grasped with the thumb and index finger and carefully pulled from the anchor elements and taken out of the mouth.

## 21.3 Cleaning and care

We recommend cleaning teeth and dentures after every meal. Cleaning of dentures includes cleaning of the connecting element. Gentlest cleaning can be achieved by cleaning the restoration under running water with a soft toothbrush and the connecting element in the mouth with an interdental brush. The most intensive cleaning of the restoration is achieved with the aid of an ultrasonic device and a cleaning additive suitable for dentures.

Never clean the high precision connecting elements with toothpaste as this could lead to damage. Caution should also be exercised in the case of aggressive cleaning agents or tablets as this could damage the high-quality connecting element or impair its function.

Regular cleaning of the anchorage can prevent inflammation of the soft tissue.

## 22 Ordering information

The information relevant to your order can be found in the product list in Section 29 of this document. The product information is also helpful. This and other relevant documents can be found on the website [www.cmsa.ch/docs](http://www.cmsa.ch/docs) by entering the relevant product name.

## 23 Availability

Some of the products described in this document may possibly not be available in all countries.

## 24 Traceability of the lot number

The lot numbers of all parts used must be documented to ensure traceability.

## 25 Complaint

Cendres+Métaux SA must be notified immediately of any incident that has occurred with regard to the product. To do this, please contact your customer advisor or send us your message by e-mail to the address [complaints-cmbrand@cmsa.ch](mailto:complaints-cmbrand@cmsa.ch). In serious cases, also send a report to the competent authority where you are domiciled.

## 26 Safe disposal

The products must be disposed of in accordance with local laws and environmental regulations, taking into account the level of contamination. Cendres+Métaux Lux SA would be very pleased to accept precious metal waste. For information and additional details, please contact your Cendres+Métaux SA representative.

## 27 Trademarks

Registered trademarks of Cendres+Métaux Holding SA, Biel/Bienne, Switzerland include:

CM LOC® / Pekkton® / Elitor® / Ceramicor®

Unless explained specifically, all products marked with "®" are not registered trademarks of Cendres+Métaux Holding SA, but registered trademarks of the respective manufacturer.

## 28 Disclaimer

The manufacturer rejects any liability for damages resulting from non-compliance with these Instructions for Use. Cendres+Métaux SA products are parts of an overall concept and may only be used or combined with the appropriate original components and instruments. Otherwise, the manufacturer rejects any responsibility and liability. In case of complaints, please always include the lot number.

The use of third party products not distributed by Cendres+Métaux SA in connection with the products mentioned in the product list in Section 29 will void any warranty or other express or implied obligation of Cendres+Métaux SA.

Responsibility regarding the suitability of a product for the specific patient case is at the discretion of the specialist.

Cendres+Métaux SA disclaims any express or implied liability and shall not be responsible for any direct, indirect, punitive or other damages arising from or in connection with errors in professional judgement or practice in the use of Cendres+Métaux SA products.

The specialist is obliged to regularly study the latest developments of the products mentioned in the product list in Section 29 and their applications.

It should be noted that the descriptions contained in this document are not sufficient for the immediate application of Cendres+Métaux SA products. Expertise in dentistry, dental technology and instructions by an experienced specialist in the use of the products mentioned in the product list under Section 29 is always necessary.

In case of inconsistencies in translations, the English language version shall prevail.

## 29 Product list

| Cat. No. | Product name  | Material               | Single use | Labelling | Basic UDI-DI      |
|----------|---|------------------------|------------|-----------|-------------------|
| 05001605 | CM LOC® Male part C, casting-on / soldering technique | Ceramicor®             | Yes        | CE 0483   | 764016651000050DW |
| 05001606 | CM LOC® Male part E, laser welding technique          | Elitor®                | Yes        | CE 0483   | 764016651000050DW |
| 05003001 | CM LOC® Basic Set Titanium                            | TiAl6 V4 ELI, (Grade5) | Yes        | CE 0483   | 764016651000057EC |
| 05001995 | CM LOC® Housing Titanium for Pekkton® Inserts         | Pekkton®               | Yes        | CE 0483   | 764016651000053E4 |
| 05001314 | CM LOC® Retention insert, extra-low                   | TiAl6 V4 ELI, (Grade5) | Yes        | CE 0483   | 764016651000053E4 |
| 05001315 | CM LOC® Retention insert, low                         | Pekkton®               | Yes        | CE 0483   | 764016651000053E4 |
| 05001316 | CM LOC® Retention insert, medium                      | Pekkton®               | Yes        | CE 0483   | 764016651000053E4 |
| 05001317 | CM LOC® Retention insert, strong                      | Pekkton®               | Yes        | CE 0483   | 764016651000053E4 |
| 05001328 | CM LOC® Processing insert                             | Pekkton®               | Yes        | CE        | 764016651000007DV |
| 05001306 | CM LOC® Housing Pekkton® for Pekkton® Inserts         | Pekkton®               | Yes        | CE 0483   | 764016651000053E4 |
| 07000201 | CM LOC® Spacer  | Pekkton®               | Yes        | CE        | 764016651000026DZ |
| 07000202 | CM LOC® Block-out spacer                              | Santoprene             | Yes        | CE        | 764016651000027E3 |
| 07000204 | CM LOC® Analog  | TiAl6 V4 ELI, (Grade5) | Yes        | CE        | 764016651000034DY |
| 07000205 | CM LOC® Multi-Tool for Pekkton® Retention insert      | TiAl6 V4 ELI, (Grade5) | No         | CE        | 764016651000001DH |
| 07000206 | CM LOC® Screw Driver                                  | TiAl6 V4 ELI, (Grade5) | No         | CE        | 764016651000022DR |
| 07000213 | CM LOC® Impression part                               | Pekkton®               | Yes        | CE        | 764016651000017DY |
| 07000217 | CM LOC® Housing Extractor                             | Steel                  | No         | CE        | 764016651000009DZ |
| 07000200 | CM LOC® Instrument set                                | n/a                    | No         | n/a       | 764016651000025DX |

## 30 Labelling on packaging/symbols



Date of manufacture



Manufacturer



Catalogue number



Lot number



Quantity

[www.cmsa.ch/docs](http://www.cmsa.ch/docs)

Observe the Instructions for Use, which are available in electronic form at the address specified.

Rx only

Attention: According to US federal law, this product may only be sold by or on behalf of a physician.



Cendres+Métaux products with CE labelling meet the requirements of the relevant European requirements.



Do not re-use



Non-sterile



Protect from sunlight



Attention, observe accompanying documents



Clear product identification



European Authorised Representative



Importer



Medical device

